




Cell Harvesting by Temperature Reduction

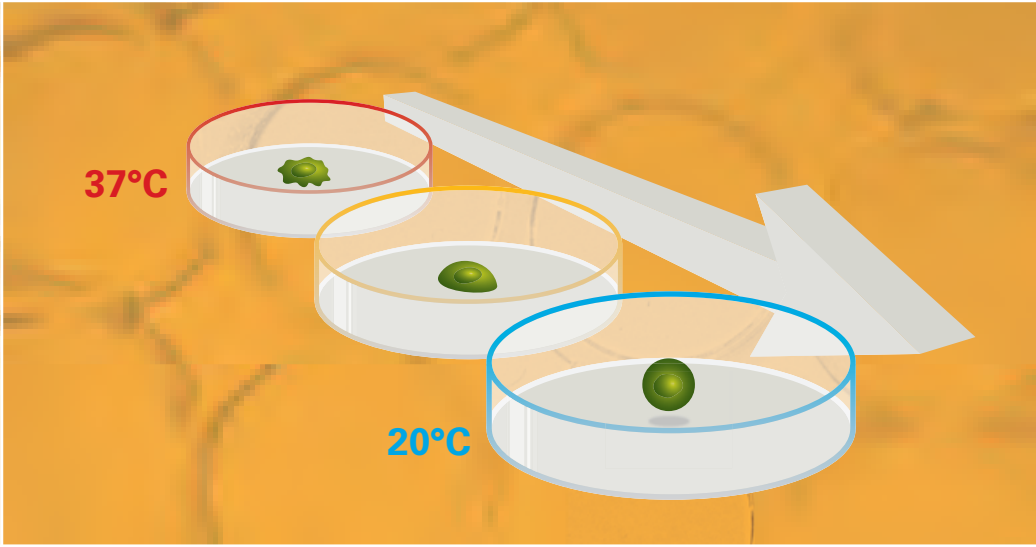

Thermo Scientific Nunc
UpCell Surface

Features and Benefits


Preserving Cell Surface Proteins

The Thermo Scientific Nunc UpCell Surface enables harvesting of cells with high viability and intact surface proteins for culture passaging, single-cell analyses and cell transplantation research.


- No trypsin – preserves cell surface proteins
- No scraping – get high cell viability
- Minimal hands-on time
- Quick, clean and simple - just reduce the temperature

Cell harvesting using enzymatic digestion, such as trypsinization, results in degradation of cell surface proteins. These proteins are important for the interactions between the cell and the environment. For example, cell surface proteins are involved in the cell's response to the extracellular matrix, to other cells and to growth factors and other soluble mediators. Some cell surface proteins are involved in the ion homeostasis of the cell, whereas other cell surface proteins are used as antigens or markers in cell analysis and enrichment procedures.

Cells harvested by trypsinization

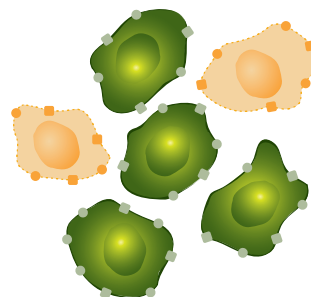

Cells harvested from UpCell Surface by temperature reduction

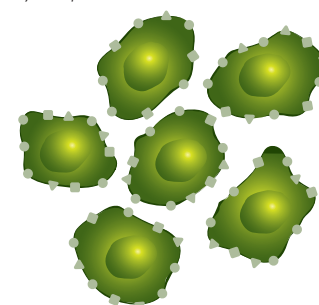


Creating 3D Tissue Models

The UpCell™ Surface enables harvesting of cell sheets and creation of tissue models held together by normal cell junctions and extracellular matrix deposited by the cells.


- No scaffold – make 3D tissue models without exogenous material
- No uneven cell distribution – control the spatial distribution of cells in 3D
- Endless possibilities for mixing cell types and creating 3D co-cultures
- Just harvest and stack cell sheets

In tissue engineering, three-dimensional (3D) tissue models or transplants are typically prepared by seeding a cell suspension on a pre-fabricated scaffold. Scaffold materials are not produced by the cells to reside in the engineered tissue, and are most often materials foreign to the body or from another species (xenogeneic), such as, poly lactic acid (PLA), poly glycolic acid (PGA), alginate, gelatin and collagen. Problems often encountered using scaffolds for tissue engineering include uneven cell distribution and difficulties in controlling the spatial distribution of different cell types. After transplantation, there can be host inflammatory reactions and fibrous tissue formation due to the exogenous scaffold material.

Tissue engineered using scaffold

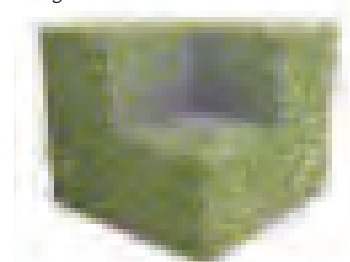


Tissue engineered using UpCell Surface

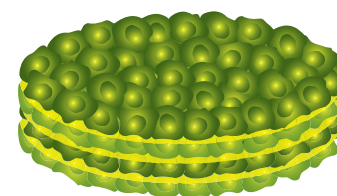

Thermo Scientific Nunc UpCell Surface

Temperature-Responsive Cell Culture Surface


UpCell Surface is designed to respond to changes in temperature. It releases adherent cells by a simple reduction of the temperature of the cell culture. Products with UpCell Surface include Thermo Scientific Nunc MicroWell plates, multidishes and dishes.

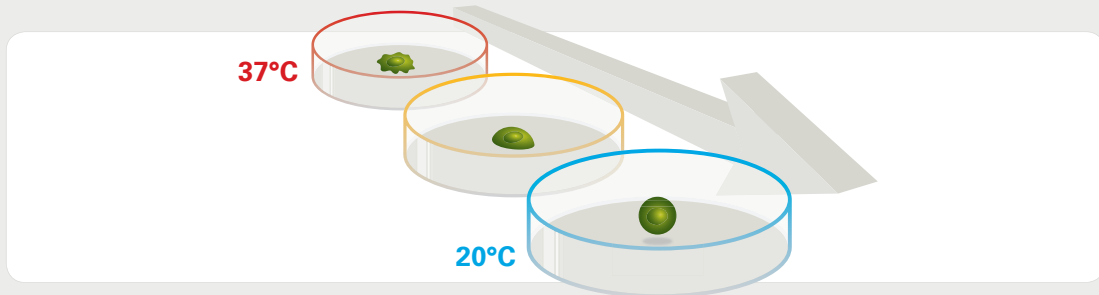

The covalently immobilized polymer poly(N-isopropylacrylamide), or PIPAAm, forms an even and thin layer on the cultureware. The PIPAAm layer is slightly hydrophobic at 37°C, allowing cells to attach and grow. When the temperature of the culture is reduced to below 32°C, the PIPAAm layer becomes very hydrophilic, binds water and swells, resulting in the release of adherent cells.

Extracellular Matrix is Harvested with the Cells

Depending on the degree of confluence of the culture, and the harvesting technique, single cells or cell sheets can be harvested from the UpCell Surface. Because the extracellular matrix under the cultured cells is harvested with the cells, cell sheets are naturally adhesive to other cell sheets and to cell surfaces in the body.

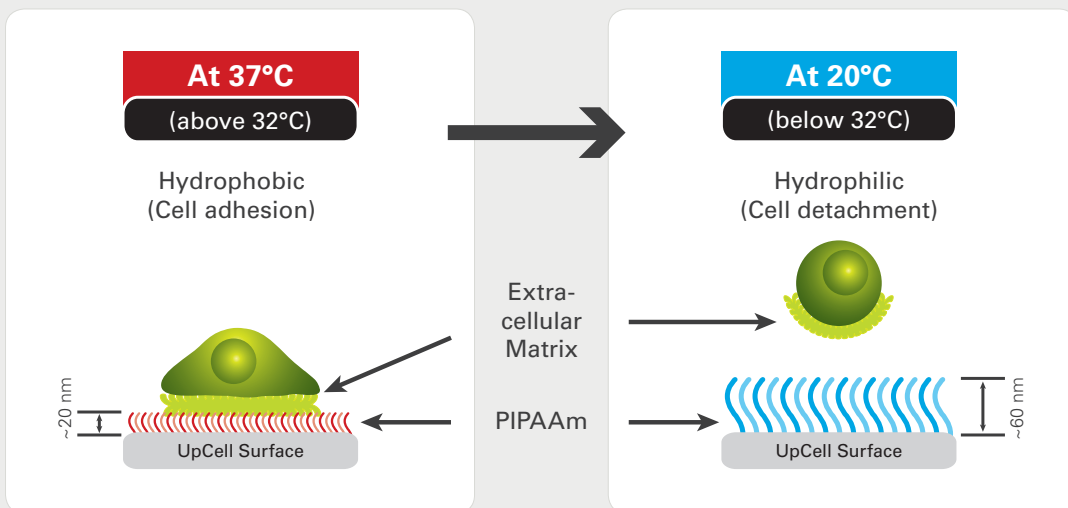

Working with Single Cells


Single-cell suspensions harvested from cultureware with Thermo Scientific Nunc UpCell Surface can be

- Analyzed *in vitro*, for example, by flow cytometry
- Manipulated to purify certain cell types
- Re-seeded in cultureware with UpCell Surface or traditional cell cultureware as part of a passaging procedure
- Used in cell transplantation research

Grow your cells in cultureware with UpCell Surface, reduce the temperature, and harvest your cells. It is that simple!

The Nunc UpCell Surface enables harvesting of cells with high viability and intact cell surface receptors and antigens. Even harvest cell types that are difficult to detach by other methods, and keep unwanted cell activation to a minimum.

Traditional cell harvesting by enzymatic and mechanical methods often compromises the integrity of surface proteins and the viability of harvested cells. By contrast, UpCell Surface allows cell harvesting by simply reducing the temperature of the cell culture, resulting in cell populations with preserved cell surface proteins and high cell viability – and there is no need for enzyme removal or inhibition.

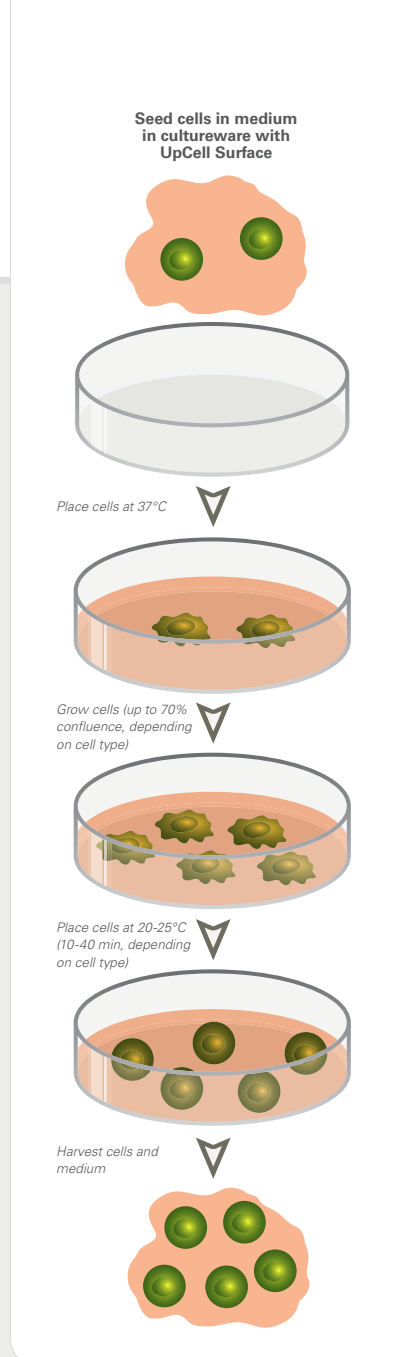

Examples of applications using single-cell suspensions harvested by temperature reduction

Cell Type	Application	Reference
Macrophages (mouse)	Detachment of cells that are otherwise difficult to detach	Application Note 1
Bone marrow cells & preadipocytes (human)	Cell surface protein preservation (flow cytometry)	Application Note 2
Microglia (rat)	Analysis (detachment and function)	Nakajima et al., 2001
Monocytes and macrophages (human)	Re-seeding/passaging	Collier et al., 2002
Monocytes and macrophages (human)	Activation & analysis (structural)	Gordon and Freedman, 2006
Basophilic cell line RBL-2H3 (rat)	Antigen-mediated degranulation measured by surface plasmon resonance	Yanase et al., 2007


Working with Cell Sheets

Cell sheets harvested from cultureware with Thermo Scientific Nunc UpCell Surface can be

- Analyzed *in vitro*, for example, by electron microscopy
- Re-plated to cultureware with UpCell Surface or traditional cell cultureware
- Used in different transplantation models
- Stacked on top of another cell sheet (see next section, Tissue Models - Cell Sheet Engineering)

The UpCell Surface enables harvesting of contiguous cell sheets with preserved cell polarization and held together by normal cell junctions and extracellular matrix. Grow your cells to confluence in cultureware with UpCell Surface, apply the supplied membrane (only Cat. Nos. 174901 6 well multidish and 174904 3.5 cm dish), reduce the temperature, and harvest the cell sheet. It is that easy!

Anchorage-dependent cells in culture initially attach to proteins adsorbed to the cultureware from the medium, and produce and deposit a subcellular matrix during the course of the cultivation (Brevig et al., 2006; Kushida et al., 1999; Pompe et al., 2003). Cells in culture are held together by cell junctions and the matrix deposited by the cells. Traditional enzymatic or mechanical harvesting disrupts these cell-to-cell and cell-to-matrix contacts, as well as the subcellular matrix, and destroys the integrity and polarization of the culture. With the UpCell Surface, cells can be detached as contiguous sheets from the cultureware. The subcellular matrix and the supplied membrane provides the mechanical strength necessary for the handling of the detached cell sheet.

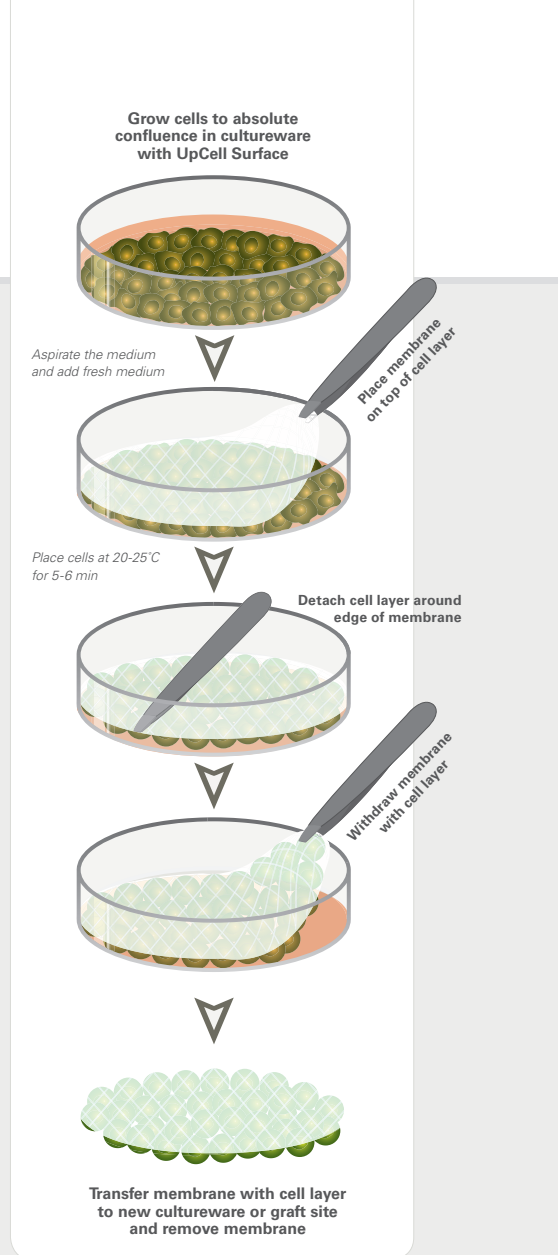

Examples of applications using cell sheets harvested by temperature reduction

Cell Type	Application	Reference
Aortic endothelial cells (cow)	Analysis (structural and matrix deposition)	Kushida et al., 1999
Keratinocytes (human)	Analysis (electron microscopy)	Yamato et al., 2001
Urothelial cells (human)	Analysis (electron microscopy)	Shiroyanagi et al., 2003
Retinal pigment epithelial cell line ARPE-19 (human)	Analysis (light microscopy)	Kubota et al., 2006
Kidney epithelial cells (human and dog)	Re-plating to traditional cultureware & analysis (electron and fluorescence microscopy)	Application Note 3 Kushida et al., 2005
Lung cells (rat)	Re-plating to traditional cultureware & analysis (fluorescence microscopy)	Nandkumar et al., 2002
Smooth muscle cells & fibroblasts (human)	Re-plating to PIPAAm surface & analysis (functional) & transplantation	Hobo et al., 2008
Mesenchymal stem cells & skin fibroblasts (rat)	Analysis (structural and functional) & transplantation	Miyahara et al., 2006
Corneal stem cells (human and rabbit)	Analysis (structural) & transplantation	Nishida et al., 2004
Corneal endothelial cells (human)	Analysis (structural and functional) & transplantation	Sumide et al., 2006
Oral mucosal epithelial cells (dog)	Analysis (structural) & transplantation	Ohki et al., 2006
Tracheal epithelial cells (rabbit)	Transplantation	Kanzaki et al., 2006
Periodontal ligament cells (human)	Transplantation	Hasegawa et al., 2005


Working with Tissue Models - Cell Sheet Engineering

Cell sheet constructs prepared in cultureware with Thermo Scientific Nunc UpCell Surface can be

- Analyzed *in vitro*, for example, by functional tissue-specific tests
- Cultivated *in vitro*, for example, as 3D co-cultures
- Used in different transplantation models, where cell sheets are stacked before or during the transplantation procedure

The UpCell Surface enables harvested cell sheets to be stacked in order to form 3D tissue models. Grow your cells to confluence in cultureware with UpCell Surface, harvest the cell sheet, and transfer the cell sheet to another cell sheet. No scaffold is needed!

Stacking of cell sheets, also known as Cell Sheet Engineering, was pioneered by Okano and colleagues (Yamada et al., 1990; Yang et al., 2005 and 2007). The preserved subcellular matrix of a harvested cell sheet provides the adhesive necessary for stacking. It functions as a natural glue to bond the cell sheet to an underlying recipient cell sheet or to a recipient site in a transplantation model, without the use of fibrin glue or sutures.

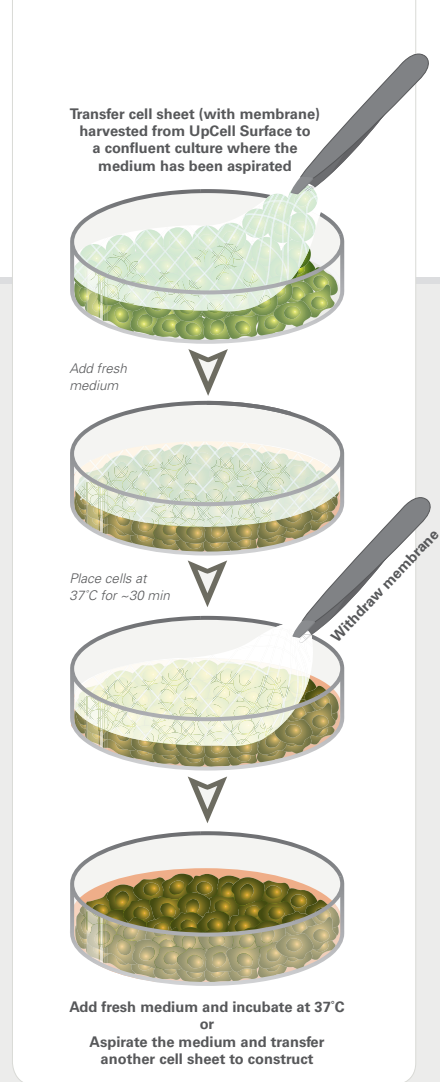

Examples of Cell Sheet Engineering using cell sheets harvested by temperature reduction

Cell Type	Application	Reference
Aortic endothelial cells (human) & hepatocytes (rat)	Cultivation (3D co-culture) & analysis (structural)	Harimoto et al., 2002
Hepatocytes (mouse and human)	Analysis (structural and functional) & stacking during transplantation	Ohashi et al., 2007
Skeletal myoblast (dog)	Analysis (structural) & stacking during transplantation	Hata et al., 2006
Lung and skin fibroblasts (rat)	Analysis (structural) & stacking during transplantation	Kanzaki et al., 2007
Cardiomyocytes (rat)	Cultivation & analysis (structural and functional) & stacking before transplantation	Sekine et al., 2006; Sekiya et al., 2006; Shimizu et al., 2002 and 2006

Quality Assurance

Quality is inherent in our culture. From product development and sourcing raw materials to manufacturing and customer service, quality is reflected in every Thermo Scientific Nunc product.

A certificate of quality is packed in each box of cultureware with Thermo Scientific Nunc UpCell Surface. This certificate is your guarantee that the product has been validated according to the following tests:

Cell growth

Each manufacturing lot is sampled and subjected to performance testing for growth with the 3T3-Swiss Albino cell line (derived from a mouse embryo fibroblast) in accordance with standard operating procedures. Acceptance level: minimum of 80% confluence.

Cell detachment

The manufacturing lot is sampled and subjected to performance testing for cell detachment by temperature reduction with the 3T3-Swiss Albino cell line in accordance with standard operating procedures. Adherent cells are detached by temperature decreasing treatment (under 32°C) and the degree of detachment is measured. Acceptance level: detachment of 50% or more of the cells.

Sterility

Sterility is obtained by using ethylene oxide gas according to ISO 11135-1 (Sterilization of health care products, Ethylene oxide, Part 1: Requirements for development, validation and routine control of a sterilization process for medical devices).

Non-pyrogenic

Representative samples of products with UpCell Surface are tested according to the principles of the LAL-test described in the FDA guidelines and certified nonpyrogenic with a documented endotoxin level of less than 20 endotoxin units/device (0.5 Endotoxin units/mL) as stated in the USP.

Toxicity

The material has successfully passed the USP biological reactivity Class VI test - 50°C (7 days implant). Cytotoxicity test according to ISO 10993-5: Biological evaluation of medical devices - Part 5: Tests for *in vitro* cytotoxicity.

Products

Cultureware with UpCell Surface

With lid. Sterile

Cat. No.	174897	174898	174899	174900	174901	174902	174905	174903	174906	174904
Format	96 MicroWell Plate w/ flat bottom	48 well Multidish	24 well Multidish	12 well Multidish	6 well Multidish	10 cm Dish	10 cm Dish with grid	6 cm Dish	6 cm Dish with grid	3.5 cm Dish
Number of wells	96	48	24	12	6	1	1	1	1	1
Culture area, cm ² /well	0.33	1.1	1.9	3.5	9.6	56.7	56.7	21.5	21.5	8.8
Max. external dimensions, mm	128x86	128x86	128x86	128x86	128x86	92x17	92x17	60x15	60x15	40x12
Suggested working volume, mL/well	0.2	0.5	1	2	3	12.5	12.5	5	5	3
Airvent	+	+	+	+	+	+	+	+	+	+
Units per pack/case	1/8	1/6	1/6	1/6	1/6	1/6	1/6	5/30	5/30	5/30

Cultureware with UpCell Surface is intended for research purposes only and single-use only. Any other use is not warranted by Thermo Fisher Scientific. Do not use the product for clinical or diagnostic purposes.

UpCell Surface is licensed from CellSeed Inc. Made in Japan.

References


- Brevig T, Røhrmann JH, Riemann H. *Oxygen reduces accumulation of type IV collagen in endothelial cell subcellular matrix via oxidative stress.* *Artif. Organs.* 2006; 30:915-921.
- Collier TO, Anderson JM, Kikuchi A, Okano T. *Adhesion behavior of monocytes, macrophages, and foreign body giant cells on poly (N-isopropylacrylamide) temperature-responsive surfaces.* *J. Biomed. Mater. Res.* 2002; 59:136-143.
- Gordon IO, Freedman RS. *Defective antitumor function of monocyte-derived macrophages from epithelial ovarian cancer patients.* *Clin. Cancer Res.* 2006; 12:1515-1524.
- Harimoto M, Yamato M, Hirose M, Takahashi C, Isoi Y, Kikuchi A, Okano T. *Novel approach for achieving double-layered cell sheets co-culture: overlaying endothelial cell sheets onto monolayer hepatocytes utilizing temperature-responsive culture dishes.* *J. Biomed. Mater. Res.* 2002; 62:464-470.
- Hasegawa M, Yamato M, Kikuchi A, Okano T, Ishikawa I. *Human periodontal ligament cell sheets can regenerate periodontal ligament tissue in an athymic rat model.* *Tissue Eng.* 2005; 11:469-478.
- Hata H, Matsumiya G, Miyagawa S, Kondoh H, Kawaguchi N, Matsuura N, Shimizu T, Okano T, Matsuda H, Sawa Y. *Grafted skeletal myoblast sheets attenuate myocardial remodeling in pacing-induced canine heart failure model.* *J. Thorac. Cardiovasc. Surg.* 2006; 132:918-924.
- Hobo K, Shimizu T, Sekine H, Shin'oka T, Okano T, Kurosawa H. *Therapeutic angiogenesis using tissue engineered human smooth muscle cell sheets.* *Arterioscler. Thromb. Vasc. Biol.* 2008; 28:637-643.
- Kanzaki M, Yamato M, Hatakeyama H, Kohno C, Yang J, Umamoto T, Kikuchi A, Okano T, Onuki T. *Tissue engineered epithelial cell sheets for the creation of a bioartificial trachea.* *Tissue Eng.* 2006; 12:1275-1283.
- Kanzaki M, Yamato M, Yang J, Sekine H, Kohno C, Takagi R, Hatakeyama H, Isaka T, Okano T, Onuki T. *Dynamic sealing of lung air leaks by the transplantation of tissue engineered cell sheets.* *Biomaterials* 2007; 28:4294-4302.
- Kubota A, Nishida K, Yamato M, Yang J, Kikuchi A, Okano T, Tano Y. *Transplantable retinal pigment epithelial cell sheets for tissue engineering.* *Biomaterials* 2006; 27:3639-3644.
- Kushida A, Yamato M, Konno C, Kikuchi A, Sakurai Y, Okano T. *Decrease in culture temperature releases monolayer endothelial cell sheets together with deposited fibronectin matrix from temperature-responsive culture surfaces.* *J. Biomed. Mater. Res.* 1999; 45:355-362.
- Kushida A, Yamato M, Isoi Y, Kikuchi A, Okano T. *A noninvasive transfer system for polarized renal tubule epithelial cell sheets using temperature-responsive culture dishes.* *Eur. Cell Mater.* 2005;10:23-30.
- Miyahara Y, Nagaya N, Kataoka M, Yanagawa B, Tanaka K, Hao H, Ishino K, Ishida H, Shimizu T, Kangawa K, Sano S, Okano T, Kitamura S, Mori H. *Monolayered mesenchymal stem cells repair scarred myocardium after myocardial infarction.* *Nat. Med.* 2006; 12:459-465.
- Nakajima K, Honda S, Nakamura Y, López-Redondo F, Kohsaka S, Yamato M, Kikuchi A, Okano T. *Intact microglia are cultured and non-invasively harvested without pathological activation using a novel cultured cell recovery method.* *Biomaterials* 2001; 22:1213-1223.
- Nandkumar MA, Yamato M, Kushida A, Konno C, Hirose M, Kikuchi A, Okano T. *Two-dimensional cell sheet manipulation of heterotypically co-cultured lung cells utilizing temperature-responsive culture dishes results in long-term maintenance of differentiated epithelial cell functions.* *Biomaterials* 2002; 23:1121-1130.
- Nishida K, Yamato M, Hayashida Y, Watanabe K, Maeda N, Watanabe H, Yamamoto K, Nagai S, Kikuchi A, Tano Y, Okano T. *Functional bioengineered corneal epithelial sheet grafts from corneal stem cells expanded ex vivo on a temperature-responsive cell culture surface.* *Transplantation* 2004; 77:379-385.
- Ohashi K, Yokoyama T, Yamato M, Kuge H, Kanehiro H, Tsutsumi M, Amanuma T, Iwata H, Yang J, Okano T, Nakajima Y. *Engineering functional two- and three-dimensional liver systems in vivo using hepatic tissue sheets.* *Nat. Med.* 2007; 13:880-885.
- Ohki T, Yamato M, Murakami D, Takagi R, Yang J, Namiki H, Okano T, Takasaki K. *Treatment of oesophageal ulcerations using endoscopic transplantation of tissue-engineered autologous oral mucosal epithelial cell sheets in a canine model.* *Gut* 2006; 55:1704-1710.
- Pompe T, Kobe F, Salchert K, Jørgensen B, Oswald J, Werner C. *Fibronectin anchorage to polymer substrates controls the initial phase of endothelial cell adhesion.* *J. Biomed. Mater. Res. A.* 2003; 67:647-657.
- Sekine H, Shimizu T, Yang J, Kobayashi E, Okano T. *Pulsatile myocardial tubes fabricated with cell sheet engineering.* *Circulation* 2006; 114 (Suppl):I87-I93.
- Sekiya S, Shimizu T, Yamato M, Kikuchi A, Okano T. *Bioengineered cardiac cell sheet grafts have intrinsic angiogenic potential.* *Biochem. Biophys. Res. Commun.* 2006; 341:573-582.
- Shimizu T, Yamato M, Isoi Y, Akutsu T, Setomaru T, Abe K, Kikuchi A, Umezumi M, Okano T. *Fabrication of pulsatile cardiac tissue grafts using a novel 3D dimensional cell sheet manipulation technique and temperature-responsive cell culture surfaces.* *Circ. Res.* 2002; 90:e40-e48.
- Shimizu T, Sekine H, Isoi Y, Yamato M, Kikuchi A, Okano T. *Long-term survival and growth of pulsatile myocardial tissue grafts engineered by the layering of cardiomyocyte sheets.* *Tissue Eng.* 2006; 12:499-507.
- Shiroyanagi Y, Yamato M, Yamazaki Y, Toma H, Okano T. *Transplantable urothelial cell sheets harvested noninvasively from temperature-responsive culture surfaces by reducing temperature.* *Tissue Eng.* 2003; 9:1005-1012.
- Sumide T, Nishida K, Yamato M, Ide T, Hayashida Y, Watanabe K, Yang J, Kohno C, Kikuchi A, Maeda N, Watanabe H, Okano T, Tano Y. *Functional human corneal endothelial cell sheets harvested from temperature-responsive culture surfaces.* *FASEB J.* 2006; 20:392-394.
- Yamada N, Okano T, Sakai H, Karikusa F, Sawasaki Y, Sakurai Y. *Thermo-responsive polymeric surfaces; control of attachment and detachment of cultured cells.* *Makromol. Chem. Rapid Commun.* 1990; 11:571-576.
- Yamato M, Utsumi M, Kushida A, Konno C, Kikuchi A, Okano T. *Thermo-responsive culture dishes allow the intact harvest of multilayered keratinocyte sheets without dispase by reducing temperature.* *Tissue Eng.* 2001; 7:473-480.
- Yanase Y, Suzuki H, Tsutsui T, Uechi I, Hiragun T, Mihara S, Hide M. *Living cell positioning on the surface of gold film for SPR analysis.* *Biosens. Bioelectron.* 2007;23:562-567.
- Yang J, Yamato M, Kohno C, Nishimoto A, Sekine H, Fukai F, Okano T. *Cell sheet engineering: recreating tissues without biodegradable scaffolds.* *Biomaterials* 2005; 26:6415-6422.
- Yang J, Yamato M, Shimizu T, Sekine H, Ohashi K, Kanzaki M, Ohki T, Nishida K, Okano T. *Reconstruction of functional tissues with cell sheet engineering.* *Biomaterials* 2007; 28:5033-5043.

Cell Culture Excellence™

Essential products for the cell culture laboratory

Our comprehensive portfolio includes advanced tools designed to help you achieve excellence at every stage of your cell culture process – from growth and passage to experimentation through characterization, analysis and storage.

To learn about our full array of cell culture products and services, go to:
www.thermo.com/cellculture

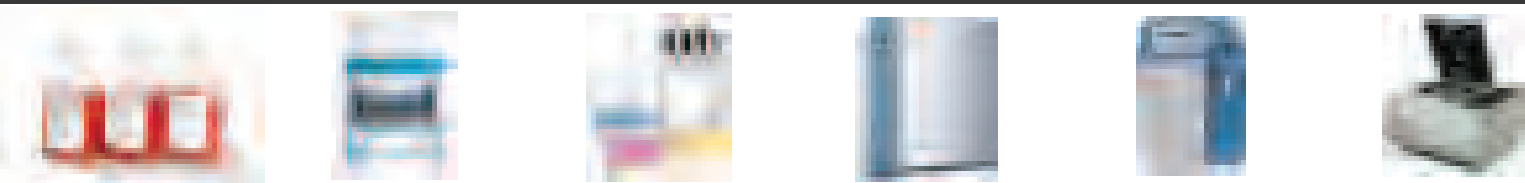

North America - Tel: 1-800-625-4327 - technical.nunc@thermofisher.com
Asia Pacific - Tel: +65 6770 2807 - intlmtg@thermofisher.com
China - Tel: 86-21-68654588 - info.nnichina@thermofisher.com
Europe (Nalgene) - Tel: +44 (0) 1432 263933 - sales@nalgene.co.uk
Europe (Nunc) - Tel: +45 4631 2000 - info.nunc@thermofisher.com
India - Tel: +91-22-67162200 - Toll Free: 1800 22 8374 - contact.lpg.in@thermofisher.com
Japan - Tel: +81 3 3816 3355 - info@nalgenunc.co.jp - www.nalgenunc.co.jp
All other locations (USA, International Department) - Tel: +1 585 899 7198 - intlmtg@thermofisher.com

www.thermo.com/upcell